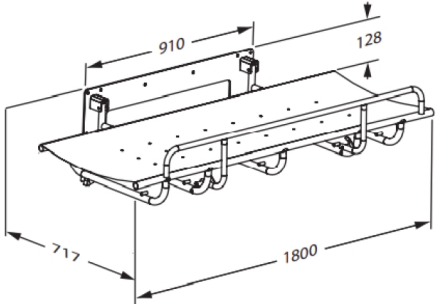




|   |  |
|---|--|
| <p><b>Product</b><br/><b>Item number</b></p>                    | <p>SCT 1000 shower change table, fixed height<br/>R8407000</p>   |
| <p><b>Colour</b><br/><br/><b>Product description</b></p>        | <p>VVS:777545380<br/>EAN:5708590294827<br/>HMI: 88217<br/>NBS:N13:40,375<br/>FDA:INT/D124320<br/>White</p>  <p>Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> <p>SCT 1000 shower change table, fixed height.<br/>Foldable. Incl. safety rail.<br/>Outer measurements 747 mm x 1800 mm.<br/>Excl. mounting kit.</p> |
| <p><b>Adjustability</b><br/><b>Product capacity</b></p>         | <p>Fixed mounting</p> <p>Maximum user weight (MUM): 150 kg.<br/>Safe working load (SWL): 150 kg.<br/>When the table is loaded with 150 kg, each mounting bolt supports 46 kg.<br/><i>This product has been tested statically with a load of 225 kg as the ISO 17966 standard requires testing with safety factor 1.5.</i></p>  |
| <p><b>Materials</b><br/><b>Ergonomy</b><br/><b>Features</b></p> | <p>The table is made from practical and durable PVC canvas which is phthalate free. The frame is powder-coated stainless steel.</p> <p>Foldable</p> <p>The integrated gas cylinders enable the table to be folded out of the way when not in use to maximise floor space in front of the unit.<br/>Splash curtain (optional) for the nursing bench protects the carer against splashing.</p>   |

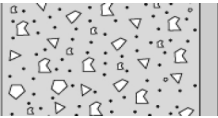
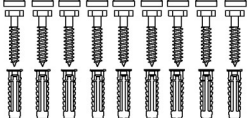

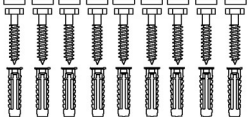
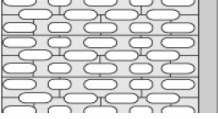
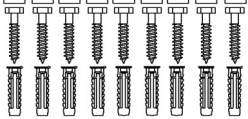
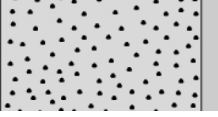
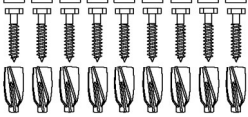
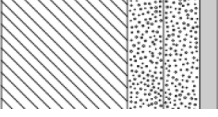
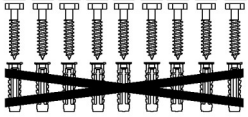


|                                  |   |
|----------------------------------|---|
| <b>Product</b>                   | SCT 1000 shower change table, fixed height  |
| <b>Item number</b>               | R8407000  |
| <b>Warranty</b>                  | Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.  |
| <b>Cleaning</b>                  | <p>SCT 1000 &amp; 2000<br/>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b><br/>Clinimax Difficil-S<br/>Erisan OXY+<br/>Sprit (husholdningsprit)<br/>Suma Bac D10<br/>Virkon<br/>Rodalon Desinfektionsspray<br/>Lysoformin Spezial<br/>CaviCide<br/>Sufanios Premium<br/>Tvättsprit 70%<br/>Accel PREvention Wipes Lingettes</p> <p><b>NON-DISINFECTANTS:</b><br/>Suma Star free D1<br/>Taski 100<br/>Taski 200<br/>Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> |
| <b>Declaration of conformity</b> | <p>CE UK MDR This product is in compliance with all relevant requirements of:<br/>- Regulation (EU) 2017/745 on medical devices.<br/><b>For the UK:</b> - Medical Device Regulation (MDR) 2002 on medical devices.</p>  |

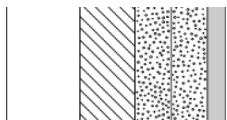
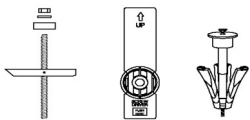
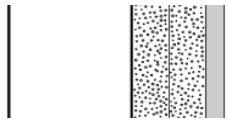
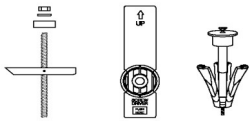
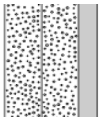



|                    |  |
|--------------------|--|
| <b>Product</b>     | SCT 1000 shower change table, fixed height |
| <b>Item number</b> | R8407000                                   |

## Installation in different wall types

|   |  |   |
|---|--|---|
|    | <b>Concrete &gt;B25</b><br>V8611 - Mounting kit, 9 hex screws Ø10 x 70 - plugs SX12.<br>To be ordered separately.  |    |
|   | <b>Solid Brick &gt;= Mz12</b><br>V8611 - Mounting kit, 9 hex screws Ø10 x 70 - plugs SX12.<br>To be ordered separately.  |   |
|  | <b>Solid Sand Lime Brick &gt;= Hlz12</b><br>V8611 - Mounting kit, 9 hex screws Ø10 x 70 - plugs SX12.<br>To be ordered separately.   |  |
|  | <b>Aerated Concrete P4</b><br>V8661 - Mounting kit, 8 bolts M5 x 50 - plugs FTP K 6.<br>To be ordered separately.  |  |
|  | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b><br>V8611 - Mounting kit, 8 screws Ø5 x 50 - plugs SX8. NB!<br>Plugs are not needed.<br>To be ordered separately. |  |



|   |   |
|---|---|
| <b>Product</b><br><b>Item number</b>  | SCT 1000 shower change table, fixed height<br>R8407000  |
|    | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b><br>Wall fasteners (9 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.<br>       |
|    | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm)</b><br>Wall fasteners (9 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.<br> |
|  | <b>Drywall (2 x 12,5 mm) water-resistant</b><br>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.<br>   |

It is the responsibility of the installer to ensure that suitable wall fasteners are used.